

ELECTROSURGICAL EQUIPMENT

P R O D U C T C A T A L O G



GENERATORS • CABLES & ADAPTERS
NON-STICK FORCEPS • ELECTRODES
SWITCHPENS • CAUTERIES

api
anthony products, inc.

Table of Contents, Company Info & Policies

Electrosurgical Generators	2-5	Electrodes	13-15
Cables, Grounding Equipment & Adapters	6	Goco Adenoid Curette & Suction.....	16
Electrosurgical Forceps	7-10	Geiger Thermal Cautery Unit.....	16
Suction Coagulators.....	11	Cauteries & Tips	17-18
Electrosurgery Pencils, Turbinate Probes.....	12	Insulated Surgical Instruments.....	19-22

Pricing Information

Please refer to our most current price list or call API at **1-800-428-1610** for the most current pricing information. You may also find prices and shop for the items shown in this catalog at www.anthonypproducts.com. We make every effort to maintain our pricing, however, due to market fluctuations, prices are subject to change without notice.

Quoting Information

We can provide quotes on products not found in this catalog. Large quantities of single products may also be quoted.

Shipping Information

- **Next Day Delivery Services:** Several options are available to meet your needs. Choose from Early Morning, Mid-morning or Afternoon*
- **2-Day Delivery Service:** Provides delivery within 2 working days.
- **Ground Delivery Service:** The most economical way to ship your package. Delivery is estimated to be within 5 working days.

NOTE: Delivery options listed are for U.S. domestic only. Large shipments requiring inside delivery and trash removal are subject to additional charges.

**Next Day Delivery times vary based on delivery provider, location, volume, etc.*

Early Morning is usually 8am-10am to most areas. Mid-morning is usually 10am-noon to most areas. Afternoon is usually 12pm-4:30pm in most areas. Residential and rural areas times may experience longer delivery times.

Warranty Information

All products in this catalog are guaranteed to be free from defects in workmanship and material. Any product which proves defective in workmanship and/or material, will either be repaired or replaced at our discretion. Warranty is void if equipment or products are not maintained properly or are not used for their intended purpose.

Minimum Order Charge

In order to continue giving the quality service you've come to expect from API, a small order handling charge of \$10.00 will be placed on orders totaling less than \$50.00.

Method of Payment

Anthony Products accepts the following forms of payment:

American Express, Discover, MasterCard and Visa.

When placing your order, please furnish the following: credit card number, CVC number, expiration date and the card holder's name.



					
Toll-Free No. 1-800-428-1610 Local No. (317) 545-6196	Fax No. (317) 543-3289	Hours of Operation Monday-Friday 8am-5pm EST	Address 7740 Records Street Indianapolis, IN 46226 USA	Online Shopping www.anthonypproducts.com E-mail api@anthonypproducts.com	

Return Goods Policy

1. A Return Merchandise Authorization (RMA#) must be obtained before returning goods. Call us for your RMA#. **Goods returned without authorization will be refused.**
2. Return Merchandise Authorization (RMA#) must be on the outside of returned package. RMA# indicates authorization to return product(s) for evaluation only. It does not imply credit, exchange, etc.
3. No goods may be returned for credit more than **10 calendar days** after the receipt of order.
4. **No credit allowed on return of special order, sterile packaged, altered or autoclaved items.**
5. All goods returned must have the original packaging. If original packaging is missing, goods are subject to a repacking fee.
6. The customer is responsible for wrapping, packing and insuring the item(s) being returned to protect from damage during shipping.
7. All returns will be **SUBJECT** to a **RESTOCK CHARGE of 25%**, unless the item(s) was shipped in error.
8. All merchandise should be insured.
9. All freight charges must be paid by the purchaser when returning goods, both to and from the purchaser.
10. **If shortages/damages occur, you have 48 hours to make a claim.** In this event, take photos to document the incident. Do not discard packing materials.
11. Orders on trial that exceed the agreement date constitutes agreement to purchase. Trial agreement must be on file.
12. All repairs sent in for evaluation/estimates must be disinfected or autoclaved before returning. Repairs that are approved for repair must be paid in advance. Please call for Evaluation/Repair Form prior to sending your product(s) to us.
13. **Return Merchandise Authorization (RMA#) expires in 10 days from issued date.**

A High Frequency Multipurpose Generator for In-Demand Physician Offices & Specialty Facilities



Bovie® Specialist|PRO Generator (Item No. ES-A1250S)

- The Bovie Specialist|PRO is the affordable solution for the busy physician office or specialty facility
- Multi-purpose unit features a modern push button design and is capable of most monopolar and bipolar surgical procedures
- Designed to operate like an operating room generator in a compact lower watt package
- Offers 5 common energy outputs controlled by the sleek, push-button design, with up to 120 watts of power.
- Modalities include Cut, Blend, Coagulation, Fulguration and Standard Bipolar
- Industry-standard connections allow for use of a variety of monopolar and bipolar accessories providing energy versatility across a variety of specialties such as Plastic Surgery, Dermatology, ER, OB/GYN, Urology, and more
- Bovie NEM™ neutral electrode monitoring increases patient safety by ensuring proper pad connectivity
- BovieDED™ (Digital Error Detection), audible alert tones and self-test circuits
- Satisfies the physician's electrosurgical needs with safety, reliability and convenience
- Uses a full line of disposable accessories for easy setup and transition between procedures

Bovie Specialist|Pro Accessories (Sold Separately)



Footswitch
(Item No. ES-BV-1253B)

Product Specifications

OUTPUT POWER

- Cut: 120 watts
- Blend: 90 watts
(60% Cut, 40% Coag.)
- Coagulation: 80 watts
- Fulguration: 40 watts
- Bipolar: 30 watts

OUTPUT FREQUENCY

- 350 - 600 kHz

LINE VOLTAGE

- 100-240 VAC

LINE FREQUENCY

- 50 - 60 Hz

WARRANTY

- 4 Years

Dimensions

- Height: 10.25" • Width 6" • Depth 12" • Weight <9lbs

Included with Purchase

- ESP1 – Sterile, Disposable, Push Button Pencil (Qty: 1)
- ES01 – Sterile, Disposable Blade Electrode (Qty: 1)
- ES02 – Sterile, Disposable Needle Electrode (Qty: 1)
- ES20 – Sterile, Disposable Ball Electrode (Qty: 1)
- ESRE – Split Disposable Return Electrode, Non Corded (Qty: 5)
- A1252C – Reusable Grounding Cable (Qty: 2)
- 110 VAC 15ft. Hospital-Grade Power Cord (Qty: 1)
- 220 VAC 15ft. Hospital-Grade Power Cord (Qty: 1)
- MC-55-237-001 – User's Guide (CD Format)

The Most Cost Effective, User Friendly Unit Available!



Bovie® Derm 942 High Frequency Desiccator (Item No. ES-A942)

- The DERM 942 and its patent pending handpiece are designed and engineered to be the most reliable and durable desiccator available today
- Unit's life is prolonged by incorporating the latest in chip technology and reducing the number of internal parts, thereby generating less internal heat
- Post-procedure and long-term cleanup are simplified by its non-absorbing, non-staining, high impact plastic case
- 40 watts of power
- Increase and decrease in 1/10th watt increments from 0-10 watts
- Digital Error Detection continuously monitors every aspect of the unit's output and instantly disables the output and displays the appropriate error code at the first sign of any problem
- The Bovie DERM 942 is a ground reference unit which means no grounding pad is required
- Less to do and less to buy makes this the most cost effective and user friendly unit on the market

Product Specifications

OUTPUT POWER

- To 40 watts (0-10 watts in 1/10th watt increments on both monopolar and bipolar)

OUTPUT FREQUENCY

- 328-492kHz

LINE VOLTAGE

- 100-240 VAC

LINE FREQUENCY

- 50 - 60 HZ

WARRANTY

- 4 Years

Dimensions

- Height: 7.4" • Width 8.98" • Depth 4.13" • Weight < 5lbs

Included with Purchase

- A902 – 3 Button Power Control Handpiece (Qty: 1 ea.)
- A804 – Non-Sterile Sharp Dermal Tips (Qty: 5 ea.)
- A805 – Sterile Sharp Dermal Tips (Qty: 2 ea.)
- A806 – Non-Sterile Blunt Dermal Tips (Qty: 5 ea.)
- A807 – Sterile Blunt Dermal Tips (Qty: 2 ea.)
- A806DE – Non-Sterile DERM-Elite™ Premium Blunt Dermal Tip (Qty: 2 ea.)
- A807DE – Sterile DERM-Elite™ Premium Blunt Dermal Tip (Qty: 2 ea.)
- A910 – Non-Sterile Disposable Handpiece Sheath (Qty: 2 ea.)
- A910ST – Sterile Disposable Handpiece Sheath (Qty: 2 ea.)
- A837 – Wall Mount Kit (Qty: 1 ea.)
- Medical-Grade Power Cords (Qty: 1 ea. 110 VAC & 220 VAC)
- MC-55-239-001 – User's Guide (CD Format)

Bovie® Derm 942 Accessories

(Sold Separately)



Footswitch
(Item No. ES-A803)



Table Top Stand
(Item No. ES-A813)



Tray for A812 Stand
(Item No. ES-A808T)



Mobile Stand
(Item No. ES-A812)

Two Machines in One Convenient Package!

Bovie® Bantam|PRO™ A952 Electrosurgical Generator (Item No. ES-A952)

- Unit is one half high frequency desiccator, an upgrade to our very popular A940, while the other half is a Bovie Generator
- All user controls and outputs on the unit's compact front panel
- Easily change power settings using the handpiece buttons or dial on the front of the unit
- Handpiece can still function even when optional footswitch is in use
- Industry standard 3/32" (0.24cm) disposable and reusable electrodes fit the A902 handpiece
- Continuously monitors the unit and instantly disables output if it detects a problem
- Post-procedure cleanup is simple (and fast) due to nonabsorbent ABS enclosures, and our use of liquid resistant membrane switches that can be wiped clean with a damp cloth



GENERATOR FEATURES

- Generator offers up to 50 watts in Cut mode, Blended Cut mode, and Pinpoint Coagulation mode
- Pad sensing technology monitors the return electrode for optimum patient protection

DESICCATOR FEATURES

- The High Frequency Desiccator is a ground referenced machine that does not require the use of a return pad
- Precision fulguration and bipolar power up to 40 watts
- Adjustable in 1/10th watt increments in the first 10 watts, for the most accurate output available

Bovie Bantam|Pro A952 Accessories

(Sold Separately)



Footswitch
(Item No. ES-A803)



Table Top Stand
(Item No. ES-A813)

Product Specifications

OUTPUT POWER

- **Cut:** 50 watts
- **Blend:** 50 watts
- **Coagulation:** 50 watts
- **Fulguration:** 40 watts
- **Bipolar:** 40 watts
- **Micro Bipolar:** 40 watts

OUTPUT FREQUENCY

- 304 - 451 kHz

LINE VOLTAGE

- 100-240 VAC

LINE FREQUENCY

- 50 Hz

WARRANTY

- 4 Years

Dimensions

- **Height:** 7.4" • **Width:** 8.98" • **Depth:** 4.13" • **Weight:** < 5lbs

Included with Purchase

- A902 – 3 Button Power Control Handpiece (Qty: 1)
- ES01 – Sterile Disposable Blade Electrodes (Qty: 3)
- ES02 – Sterile Disposable Needle Electrode (Qty: 1)
- ES20 – Sterile Disposable Ball Electrode (Qty: 1)
- A804 – Non-Sterile Sharp Dermal Tips (Qty: 5)
- A805 – Sterile Sharp Dermal Tips (Qty: 2)
- A806 – Non-Sterile Blunt Dermal Tips (Qty: 5)
- A807 – Sterile Blunt Dermal Tips (Qty: 2)
- A806DE – Non-Sterile DERM-Elite™ Premium Blunt Dermal Tip (Qty: 2)
- A807DE – Sterile DERM-Elite™ Premium Blunt Dermal Tip (Qty: 2)
- A910 – Non-Sterile Disposable Handpiece Sheath (Qty: 2)
- A910ST – Sterile Disposable Handpiece Sheath (Qty: 2)
- A837 – Wall Mount Kit (Qty: 1)
- A1252C – Reusable Grounding Cords (Qty: 2)
- ESRE – Disposable Split Grounding Pads (Qty: 5)
- Medical-Grade Power Cord 110 VAC (Qty: 1)
- MC-55-239-001 – User's Guide (CD Format)

An Economical High Frequency Desiccator for In-Office Minor Skin Procedures



Bovie® Derm 102 High Frequency Desiccator (Item No. ES-DERM-102)

- Offers ten watts of monopolar or bipolar *power, adjustable in 1/10th watt increments to desiccate and fulgurate tissue
- Designed to give practitioners the ability to perform in-office, minor skin procedures such as the removal of skin tags, benign lesions and premalignant lesions
- Power can be adjusted using the 3-button power control handpiece (A902) or using the rotary dial on the unit
- Features a large, bright, LED display for excellent visibility across the procedure room
- Comes with a total of 18 disposable dermal tips
- Affordably priced to minimize facility's expense and maximize return on investment

**To utilize this bipolar feature the A827V bipolar cable, a bipolar forcep and footswitch activation are required*

Product Specifications

OUTPUT POWER

- To 10 watts (0-10 watts in 1/10th watt increments on both monopolar and bipolar)

OUTPUT FREQUENCY

- 550 kHz

LINE VOLTAGE

- 100-240 VAC

LINE FREQUENCY

- 50 - 60 Hz

WARRANTY

- 4 Years

Dimensions

- Height: 8.2" • Width 7.2" • Depth 2.7" • Weight <3lbs

Included with Purchase

- A902 – 3 Button Power Control Handpiece (Qty: 1 ea.)
- A804 – Non-Sterile Sharp Dermal Tips (Qty: 5 ea.)
- A805 – Sterile Sharp Dermal Tips (Qty: 2 ea.)
- A806 – Non-Sterile Blunt Dermal Tips (Qty: 5 ea.)
- A807 – Sterile Blunt Dermal Tips (Qty: 2 ea.)
- A806DE – Non-Sterile DERM-Elite™ Premium Blunt Dermal Tip (Qty: 2 ea.)
- A807DE – Sterile DERM-Elite™ Premium Blunt Dermal Tip (Qty: 2 ea.)
- A910 – Non-Sterile Disposable Handpiece Sheath (Qty: 2 ea.)
- A837 – Wall Mount Kit (Qty: 1 ea.)
- Medical-Grade Power Cord 110 VAC (Qty: 1)
- MC-55-173-001 – User's Guide (CD Format)

Also Included with Purchase

- ES-A803 – Footswitch (Qty: 1)
- ES0991E – Grounding Plate (Qty: 1)
- ES8009 – Return Cord (Qty: 1)
- **Bipolar Forceps of your Choice:** see pages 7-10 (Qty: 1)
- ES-870 – Bipolar Cable (Qty: 1)
- ES1725– Needle Electrode (Qty: 1)
- ES1735– 3/4" Scalpel Electrode (Qty: 1)
- ES1764– Ball Electrode (Qty: 1)
- ES1741– 3/8" Loop Electrode (Qty: 1)

Electrosurgical Cables



Monopolar Cable
(Item No. ES470E)
Length: 10ft.



Bipolar Cable
(Item No. ES870E)
Length: 10ft.



European Bipolar Cable
(Item No. ES872E)
Length: 10ft.



Ellman Bipolar Cable
(Item No. ES19-190)
Length: 10ft.



Handswitching Cable
(Item No. ES870E-HS)
Length: 10ft.



Coated Frazier Suction Cable
(Item No. ESFS-C)
Length: 10ft.

Grounding Equipment



Grounding Plate & Cable
(Item No. ES-A802EU)
Compatible with DERM 101, DERM 102 & A942



Grounding Wristlet
(Item No. ES1391E)
Compatible with ES8008E



Return Cord
(Item No. ES8008E)
Compatible with ES1391E

Electrosurgical Terminology

Electrosurgery utilizes the heating effect of high-frequency current passing through tissue to fulgurate, desiccate, coagulate or cut.

Below are some brief descriptions of these terms.

- **Fulguration** – A technique used for superficial dehydration or coagulation of tissue. The electrode is held a short distance away and sparks to jump to the tissue.
- **Dessication** – A technique used for superficial dehydration and deliberate destruction of tissue. The electrode is placed in contact with tissue and left there to dehydrate the tissue.
- **Coagulation** – This is the sealing of small blood vessels. The electrode is left in contact with tissue for a period of time until a deep white coagulum is formed.
- **Cutting** – The rapid volatilization of cell tissue caused by a high intensity arc generated between the electrode and tissue.

Standard “Banana Plug”

(Item No. ES960A)




Connects disposable foot activated instruments to the following:

- Acmitron
- Birtcher Blendtome
- Bovie “Specialist”
- Ellman
- Malis
- Martin
- Wolf
- Most European ES units

Required for Suction Cauteries & Disposable Pens

Electrosurgical Forceps

Introducing API's new line of Non-Stick Electrosurgical Forceps...

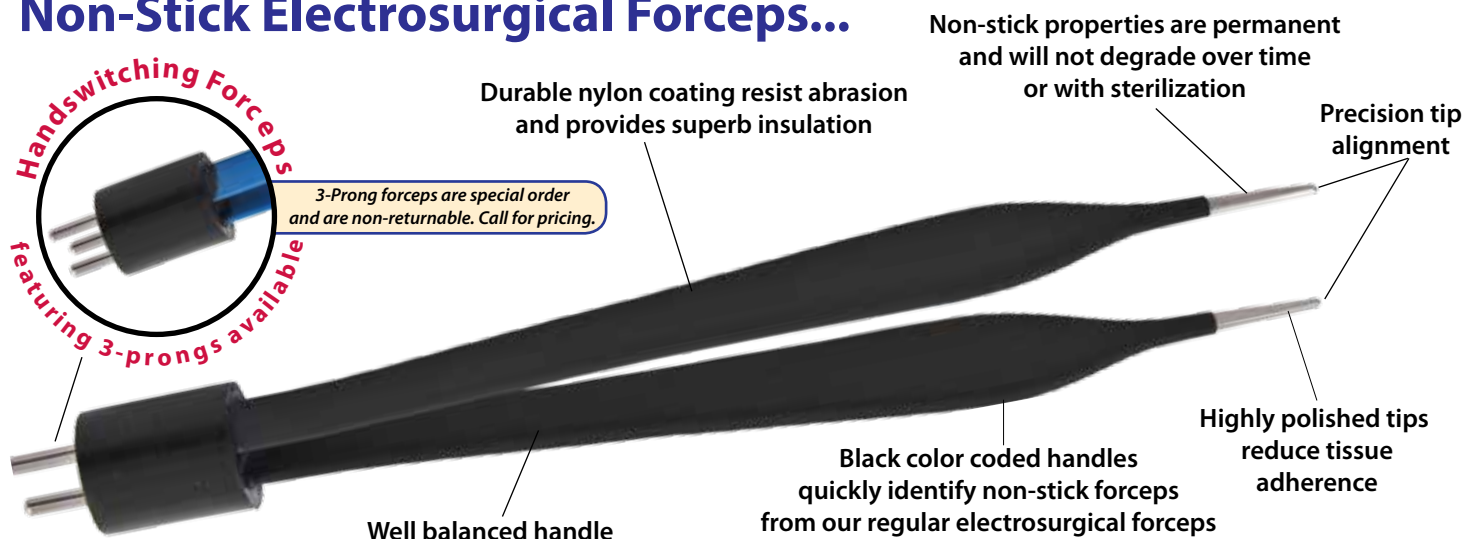

Electrosurgical Cable Options

Each electrosurgical forceps purchase includes one cable of your choice at no charge!



ES470E Monopolar Cable



ES870E Bipolar Cable



ES870E-HS Handswitching Cable



ES872E European Bipolar



ES19-190 Bipolar Cable

4 3/4" Adson Electrosurgical Forceps

Insulated, 1.0mm tip, straight



4 3/4", Non-Stick Adson Straight Black Forceps



4 3/4", Regular Adson Straight Blue Forceps

Item No.	Description	Item No.	Description
ES-200-NS	Bipolar Non-Stick Adson Forceps	ES-200	Bipolar Regular Adson Forceps
ES-200M-NS	Monopolar Non-Stick Adson Forceps	ES-200M	Monopolar Regular Adson Forceps
ES-200-BS-NS	Handswitching (Bipolar) Non-Stick Adson Forceps*	ES-200-BS	Handswitching (Bipolar) Regular Adson Forceps*
ES-200-MS-NS	Handswitching (Monopolar) Non-Stick Adson Forceps	ES-200-MS	Handswitching (Monopolar) Regular Adson Forceps



ES-200L-ANG-NS

6" Adson Electrosurgical Forceps

Insulated, 1.0mm tip



6", Non-Stick Adson Black Forceps



6", Regular Adson Straight Blue Forceps

Item No.	Description	Item No.	Description
ES-200L-NS	Bipolar Non-Stick Adson Forceps, straight tip	ES-200L	Bipolar Regular Adson Forceps, straight tip
ES-200L-ANG-NS	Bipolar Non-Stick Adson Forceps, angled tips		

#5 Jeweler Electrosurgical Forceps

4 1/2", insulated, straight



Non-Stick #5 Jeweler Straight Black Forceps

Regular #5 Jeweler Straight Blue Forceps

Item No.	Description
ES-140-NS	Bipolar Non-Stick #5 Jeweler Forceps
ES-140M-NS	Monopolar Non-Stick #5 Jeweler Forceps
ES-140-BS-NS	Handswitching (Bipolar) Non-Stick #5 Jeweler Forceps*
ES-140-MS-NS	Handswitching (Monopolar) Non-Stick #5 Jeweler Forceps

Item No.	Description
ES-140	Bipolar Regular #5 Jeweler Forceps
ES-140M	Monopolar Regular #5 Jeweler Forceps
ES-140-BS	Handswitching (Bipolar) Regular #5 Jeweler Forceps*
ES-140-MS	Handswitching (Monopolar) Regular #5 Jeweler Forceps

#7 Jeweler Electrosurgical Forceps

4 1/2", insulated, curved



Non-Stick #7 Jeweler Curved Black Forceps

Regular #7 Jeweler Curved Blue Forceps

Item No.	Description
ES-170-NS	Bipolar Non-Stick #7 Jeweler Forceps
ES-170M-NS	Monopolar Non-Stick #7 Jeweler Forceps
ES-170-BS-NS	Handswitching (Bipolar) Non-Stick #7 Jeweler Forceps*
ES-170-MS-NS	Handswitching (Monopolar) Non-Stick #7 Jeweler Forceps

Item No.	Description
ES-170	Bipolar Regular #7 Jeweler Forceps
ES-170M	Monopolar Regular #7 Jeweler Forceps
ES-170-BS	Handswitching (Bipolar) Regular #7 Jeweler Forceps*
ES-170-MS	Handswitching (Monopolar) Regular #7 Jeweler Forceps

Iris Straight Electrosurgical Forceps

4 1/2", insulated, 1.0mm tip, straight



Non-Stick Iris Straight Black Forceps

Regular Iris Straight Blue Forceps

Item No.	Description
ES-180-NS	Bipolar Non-Stick Iris Straight Forceps
ES-180M-NS	Monopolar Non-Stick Iris Straight Forceps
ES-180-BS-NS	Handswitching (Bipolar) Non-Stick Iris Straight Forceps*
ES-180-MS-NS	Handswitching (Monopolar) Non-Stick Iris Straight Forceps

Item No.	Description
ES-180	Bipolar Regular Iris Straight Forceps
ES-180M	Monopolar Regular Iris Straight Forceps
ES-180-BS	Handswitching (Bipolar) Regular Iris Straight Forceps*
ES-180-MS	Handswitching (Monopolar) Regular Iris Straight Forceps

Iris Curved Electrosurgical Forceps

4 1/2", insulated, 1.0mm tip, curved



Non-Stick Iris Curved Black Forceps

Regular Iris Curved Blue Forceps

Item No.	Description
ES-190-NS	Bipolar Non-Stick Iris Curved Forceps
ES-190M-NS	Monopolar Non-Stick Iris Curved Forceps
ES-190-BS-NS	Handswitching (Bipolar) Non-Stick Iris Curved Forceps*
ES-190-MS-NS	Handswitching (Monopolar) Non-Stick Iris Curved Forceps

Item No.	Description
ES-190	Bipolar Regular Iris Curved Forceps
ES-190M	Monopolar Regular Iris Curved Forceps
ES-190-BS	Handswitching (Bipolar) Regular Iris Curved Forceps*
ES-190-MS	Handswitching (Monopolar) Regular Iris Curved Forceps

Electrosurgical Forceps

Lucae Bayonet Electrosurgical Forceps

5 1/2", insulated, straight

Available Tip Sizes: 0.5mm, 1.0mm (stock), 1.5mm



Non-Stick Lucae Straight Black Forceps

Regular Lucae Straight Blue Forceps

Item No.	Description	Item No.	Description
ES-340-NS	Bipolar Non-Stick Lucae Forceps	ES-340	Bipolar Regular Lucae Forceps
ES-340M-NS	Monopolar Non-Stick Lucae Forceps	ES-340M	Monopolar Regular Lucae Forceps
ES-340-BS-NS	Handswitching (Bipolar) Non-Stick Lucae Forceps*	ES-340-BS	Handswitching (Bipolar) Regular Lucae Forceps*
ES-340-MS-NS	Handswitching (Monopolar) Non-Stick Lucae Forceps	ES-340-MS	Handswitching (Monopolar) Regular Lucae Forceps

Jansen Bayonet Electrosurgical Forceps

Insulated, straight

Available Tip Sizes: 0.5mm, 1.0mm (stock), 1.5mm



6" Non-Stick Jansen Bayonet Straight Black Forceps

6" Regular Jansen Bayonet Straight Blue Forceps

Item No.	Description	Item No.	Description
ES-350-NS	Bipolar Non-Stick Jansen Bayonet Forceps	ES-350	Bipolar Regular Jansen Bayonet Forceps
ES-350M-NS	Monopolar Non-Stick Jansen Bayonet Forceps	ES-350M	Monopolar Regular Jansen Bayonet Forceps
ES-350-BS-NS	Handswitching (Bipolar) Non-Stick Jansen Bayonet Forceps*	ES-350BS	Handswitching (Bipolar) Regular Jansen Bayonet Forceps*
ES-350-MS-NS	Handswitching (Monopolar) Non-Stick Jansen Bayonet Forceps	ES-350-MS	Handswitching (Monopolar) Regular Jansen Bayonet Forceps

7" Non-Stick Jansen Bayonet Straight Black Forceps

7" Regular Jansen Bayonet Straight Blue Forceps

Item No.	Description	Item No.	Description
ES-360-NS	Bipolar Non-Stick Jansen Bayonet Forceps	ES-360	Bipolar Regular Jansen Bayonet Forceps
ES-360M-NS	Monopolar Non-Stick Jansen Bayonet Forceps	ES-360M	Monopolar Regular Jansen Bayonet Forceps
ES-360-BS-NS	Handswitching (Bipolar) Non-Stick Jansen Bayonet Forceps*	ES-360BS	Handswitching (Bipolar) Regular Jansen Bayonet Forceps*
ES-360-MS-NS	Handswitching (Monopolar) Non-Stick Jansen Bayonet Forceps	ES-360-MS	Handswitching (Monopolar) Regular Jansen Bayonet Forceps

8" Non-Stick Jansen Bayonet Straight Black Forceps

8" Regular Jansen Bayonet Straight Blue Forceps

Item No.	Description	Item No.	Description
ES-370-NS	Bipolar Non-Stick Jansen Bayonet Forceps	ES-370	Bipolar Regular Jansen Bayonet Forceps
ES-370M-NS	Monopolar Non-Stick Jansen Bayonet Forceps	ES-370M	Monopolar Regular Jansen Bayonet Forceps
ES-370-BS-NS	Handswitching (Bipolar) Non-Stick Jansen Bayonet Forceps*	ES-370BS	Handswitching (Bipolar) Regular Jansen Bayonet Forceps*
ES-370-MS-NS	Handswitching (Monopolar) Non-Stick Jansen Bayonet Forceps	ES-370-MS	Handswitching (Monopolar) Regular Jansen Bayonet Forceps

Hardy Bayonet Electrosurgical Forceps

8 1/2", insulated, 1.0mm tip, angled



Non-Stick Hardy Bayonet Black Forceps

Regular Hardy Bayonet Blue Forceps

Item No.	Description	Item No.	Description
ES-470-NS	Bipolar Non-Stick Hardy Bayonet Forceps	ES-470	Bipolar Regular Hardy Bayonet Forceps
ES-470M-NS	Monopolar Non-Stick Hardy Bayonet Forceps	ES-470M	Monopolar Regular Hardy Bayonet Forceps
ES-470-BS-NS	Handswitching (Bipolar) Non-Stick Hardy Bayonet Forceps*	ES-470-BS	Handswitching (Bipolar) Regular Hardy Bayonet Forceps*
ES-470-MS-NS	Handswitching (Monopolar) Non-Stick Hardy Bayonet Forceps	ES-470-MS	Handswitching (Monopolar) Regular Hardy Bayonet Forceps

Gerald Straight Electrosurgical Forceps

7", insulated, straight

Available Tip Sizes: 1.0mm (stock), 1.5mm



Non-Stick Gerald Straight Black Forceps

Regular Gerald Straight Blue Forceps

Item No.	Description
ES-300-NS	Bipolar Non-Stick Gerald Straight Forceps
ES-300M-NS	Monopolar Non-Stick Gerald Straight Forceps
ES-300-BS-NS	Handswitching (Bipolar) Non-Stick Gerald Straight Forceps*
ES-300-MS-NS	Handswitching (Monopolar) Non-Stick Gerald Straight Forceps

Item No.	Description
ES-300	Bipolar Regular Gerald Straight Forceps
ES-300M	Monopolar Regular Gerald Straight Forceps
ES-300-BS	Handswitching (Bipolar) Regular Gerald Straight Forceps*
ES-300-MS	Handswitching (Monopolar) Regular Gerald Straight Forceps

Gerald Curved Electrosurgical Forceps

7", insulated, curved

Available Tip Sizes: 1.0mm (stock), 1.5mm



Non-Stick Gerald Curved Black Forceps

Regular Gerald Curved Blue Forceps

Item No.	Description
ES-300C-NS	Bipolar Non-Stick Gerald Curved Forceps
ES-300CM-NS	Monopolar Non-Stick Gerald Curved Forceps
ES-300C-BS-NS	Handswitching (Bipolar) Non-Stick Gerald Curved Forceps*
ES-300C-MS-NS	Handswitching (Monopolar) Non-Stick Gerald Curved Forceps

Item No.	Description
ES-300C	Bipolar Regular Gerald Curved Forceps
ES-300CM	Monopolar Regular Gerald Curved Forceps
ES-300C-BS	Handswitching (Bipolar) Regular Gerald Curved Forceps*
ES-300C-MS	Handswitching (Monopolar) Regular Gerald Curved Forceps

Cushing Electrosurgical Forceps

6 1/4", insulated, straight

Available Tip Sizes: 1.0mm (stock), 1.5mm



Non-Stick Cushing Straight Black Forceps

Regular Cushing Straight Blue Forceps

Item No.	Description
ES-230-NS	Bipolar Non-Stick Cushing Forceps
ES-230M-NS	Monopolar Non-Stick Cushing Forceps
ES-230-BS-NS	Handswitching (Bipolar) Non-Stick Cushing Forceps*
ES-230-MS-NS	Handswitching (Monopolar) Non-Stick Cushing Forceps

Item No.	Description
ES-230	Bipolar Regular Non-Stick Cushing Forceps
ES-230M	Monopolar Regular Non-Stick Cushing Forceps
ES-230-BS	Handswitching (Bipolar) Regular Non-Stick Cushing Forceps*
ES-230-MS	Handswitching (Monopolar) Regular Non-Stick Cushing Forceps

Semkin Electrosurgical Forceps

5 1/2", insulated, 0.5mm tip, straight



Item No.	Description
ES-220	Bipolar Regular Semkin Forceps
ES-220M	Monopolar Regular Semkin Forceps

Semkin Forceps not available in Non-Stick

#4 Jeweler Electrosurgical Forceps

4 1/2", insulated, 1.0mm tip, straight




Item No.	Description
ES-150	Bipolar Regular #4 Jeweler Forceps
ES-150M	Monopolar Regular #4 Jeweler Forceps
ES-150-MS	Handswitching (Monopolar) Regular #4 Jeweler Forceps

#4 Jeweler Forceps not available in Non-Stick

* Bipolar Handswitching fits Valley Lab, Conmed and Kirwan only.

Suction Coagulators, Pencil Suction Attachment

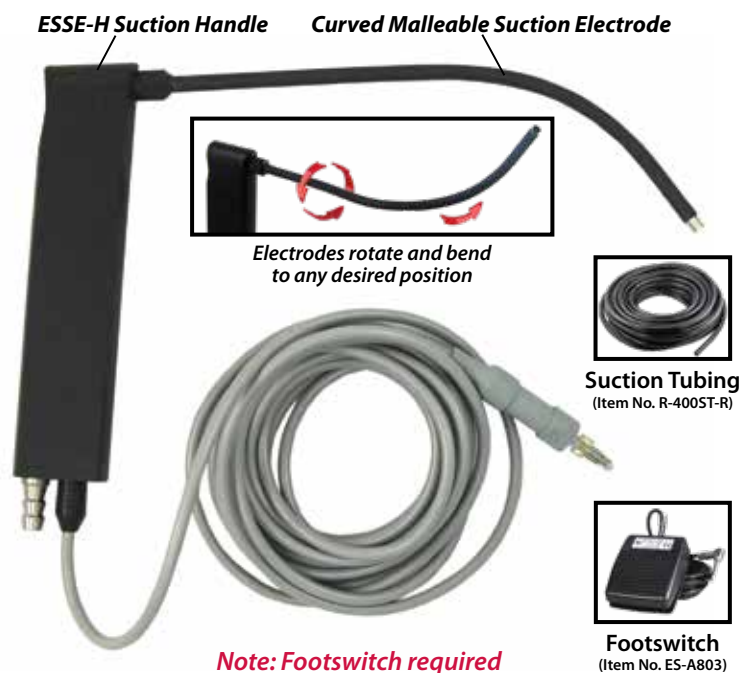

Reusable Malleable Tip Suction Electrodes

- Coagulates surgical area while simultaneously removing harmful materials including smoke, blood and mucous
- Suction coagulation handle features molded finger grips and has a suction tube and monopolar cord connection on bottom
- Can be autoclaved

Item No. Description

ESSE-H Suction Handle

ES8008E Monopolar Cord (8ft.)

R-400ST-R Suction Tubing (10ft.)

MALLEABLE SUCTION ELECTRODE TIP CHOICES

Item No. Tip Diameter Length

ESSE-8-8M 8 French (2.67mm) 8"

ESSE-10-8M 10 French (3.33mm) 8"

ESSE-12-8M 12 French (4mm) 8"

Other tips sizes available. Please call for complete list.

Bovie Suction Coagulators

- Compatible with most generator models
- 6" flexible shaft
- 8, 10 or 12 French diameters
- Choose from foot or hand activated suction coagulators
- Allows suction with or without simultaneous coagulation
- Single use/individually packaged sterile

HANDSWITCHING

Item No. Tip Diameter Length

ES-SCH08-1 8 French (2.67mm) 3m

ES-SCH10-1 10 French (3.33mm) 3m

ES-SCH12-1 12 French (4mm) 3m



FOOTSWITCHING

Item No. Tip Diameter Length

ES-SCF08-1 8 French (2.67mm) 3m

ES-SCF10-1 10 French (3.33mm) 3m

ES-SCF12-1 12 French (4mm) 3m

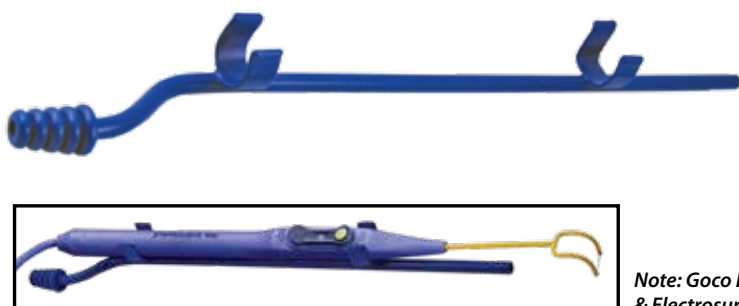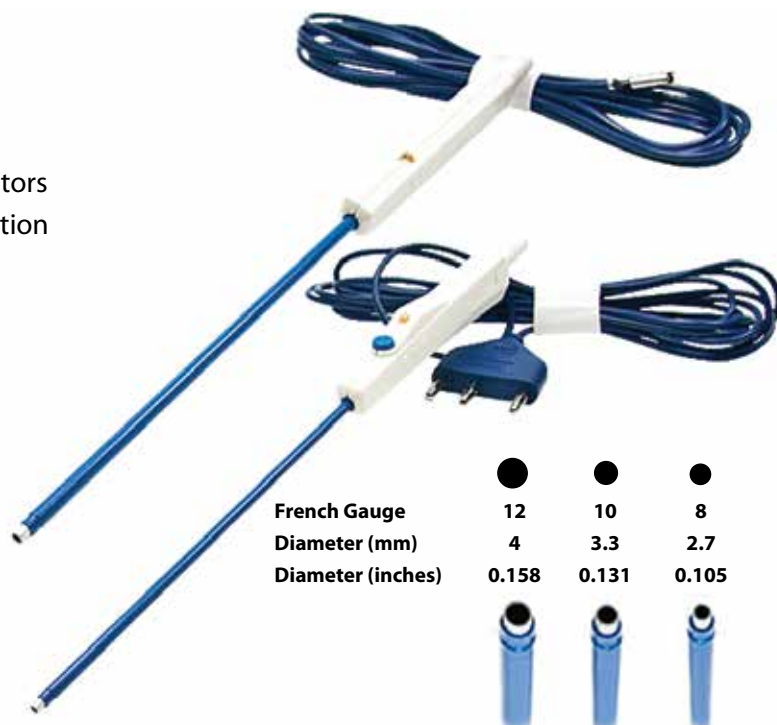

Coated Suction Attachment for Electrosurgery Pencil

(Item No. ES-FORN-LSA-LB)

- Attaches to most electrosurgery pencils
- Lightweight, ergonomically design
- Autoclavable

Electrosurgery Pencils, Bipolar Turbinate Probes

Reusable Electrosurgery Pencil

(Item No. ES-ESPR)

- Validated to 100 autoclave cycles
- Button switch activation
- Standard 3-prong connector with 10ft. attached cable
- Accepts standard 3/32" shaft electrodes (sold separately)
- Supplied non-sterile
- 1 per box

Compatible Generators

Aspen Labs, Birtcher, Bovie, Cameron Miller, Concept, Conmed, EMS/Bard, Ellman, Erbe, Eschmann, Martin, 3M-600, Neomed and Valleylab and any other generators that utilize the standard 3-terminal receiver.



Note: Hand activated generator required
(sold separately)


Disposable Electrosurgery Pencil

(Item No. ES-ESP6)

- Rocker switch activation
- Standard 3-prong connector with 10ft. attached cable
- Accepts standard 3/32" shaft electrodes
- One blade electrode included with purchase
- Packaged individually, sterile

Compatible Generators

Aspen Labs, Birtcher, Bovie, Cameron Miller, Concept, Conmed, EMS/Bard, Ellman, Erbe, Eschmann, Martin, 3M-600, Neomed and Valleylab and any other generators that utilize the standard 3-terminal receiver.

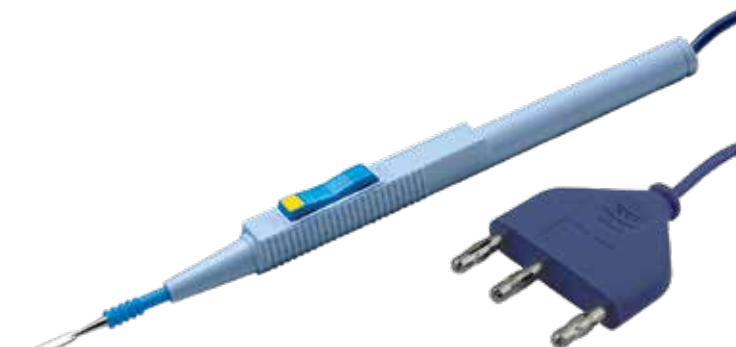

Note: Hand activated generator required
(sold separately)

Reusable Monopolar Cord & Handle

(Item No. ES1786E)

- A long-life reusable monopolar handle with 9ft. attached cable
- Accepts standard 3/32" shaft electrodes
- One blade electrode included with purchase
- Supplied non-sterile
- 1 per package



Note: Footswitch required
(sold separately)

Reusable Bipolar Turbinate Probes

- Designed for bipolar submucosal coagulation
- Two 5cm needles spaced 3mm apart
- Each probe includes one ES870E bipolar cable

Item No.	Description
ES1799E	Bayonet Bipolar Turbinate Probe
ES1799E-ANG	45° Angled Bipolar Turbinate Probe



Includes one
ES870E Bipolar Cable





Reusable Blade Electrodes

Item No.	Blade	Length
(1) ES1735	20mm	2 3/4"
(2) ES660-4	20mm	4"
(3) ES660-5	20mm	5"
(4) ES660-7	20mm	7"
(5) ES660-8	20mm	8"



Reusable Needle Electrodes

Item No.	Length
(6) ES661S	1" (standard)
(7) ES661B	1" (blunt)
(8) ES661F	1" (fine)
(9) ES661-4	4" (standard)
(10) ES661-6	6" (standard)
(11) ES661-7	7" (standard)

Reusable Electrode Features

- Standard 3/32" (2.3 mm) stainless steel shaft
- Molded safety grip insulators
- Easy insertion/removal from the surgical pencil
- Individually packaged
- Non-sterile




Reusable Ball Electrodes

Item No.	Diameter	Length
(12) ES1764	1/4"	1"
(13) ES662-4	3/16"	4"

Reusable Loop Electrodes

Item No.	Diameter
(14) ES663-1	1/4"
(15) ES663-2	5/16"
(16) ES663-3	3/8"
(17) ES663-4	1/2"

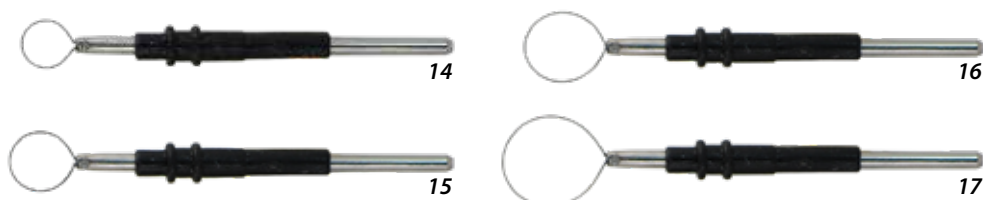

Reusable Extra Fine Wire Electrode

Item No.	Description
(18) ES666A	Angled



Reusable Teflon Coated Tungsten Fine Needle Electrodes

Item No.	Size
(19) ES04-081	5cm
(20) ES04-082	7cm
(21) ES04-083	9cm
(22) ES04-084	10cm
(23) ES04-085	13cm
(24) ES04-086	15cm



Reusable Electrode Extensions

Item No.	Size
(25) ES675-4	4"
(26) ES675-5	5"
(27) ES675-6	6"



Monopolar Brow Lift Cautery

(Item No. ES17-150)



- 7" long malleable shaft with sharp hockey stick hook
- Compatible with ES1786E Monopolar Cord with attached handle

Reusable Economical Teflon Coated Indiana Style Needles


Item No.	Description
(29) ES664	1 1/16", short, angled
(30) ES664L	3", long, angled turbinate
(31) ES664S	1 1/16", short, straight
(32) ES665	2", x-fine, straight
(33) ES665A	2 3/16", x-fine, angled



Non-Stick Reusable Electrodes

Non-Stick Reusable Electrode Features

- Standard 3/32" (2.3 mm) stainless steel shaft
- Molded safety grip insulators
- Easy insertion/removal from the surgical pencil
- Individually packaged
- Non-sterile

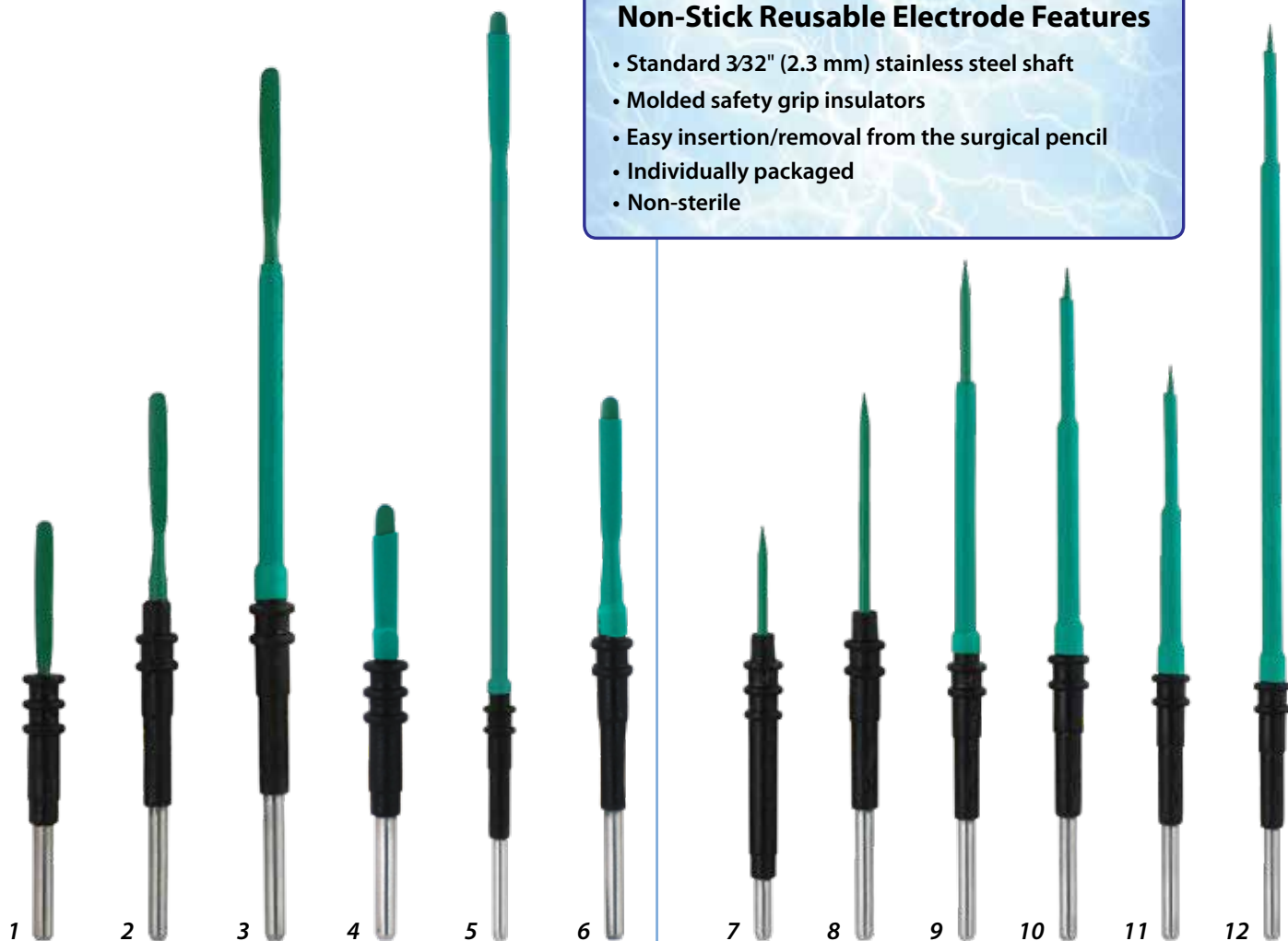

Non-Stick Blade Electrodes

	Item No.	Blade	Length
(1)	ES04-051	20mm	5cm
(2)	ES04-052	20mm	7cm
(3)	ES04-053	20mm	10cm
(4)	ES04-061	3mm	5cm
(5)	ES04-061-15	3mm	15cm
(6)	ES04-062	3mm	7cm

Non-Stick Needle Electrodes

	Item No.	Length
(7)	ES04-065	5cm
(8)	ES04-055	7cm
(9)	ES04-056	10cm
(10)	ES04-066	10cm
(11)	ES04-066-8	8cm
(12)	ES04-066-15	15cm



Non-Stick Ball Electrodes

	Item No.	Length
(13)	ES04-068	15cm
(14)	ES04-058	5cm

Adenoid Curette, Suction Catheter, Thermal Cautery Unit

Goco Electrocautery Adenoid Curette


Why cauterize when you can cut and cauterize using the GEAC tip!



- Decreased adenoid bleeding
- Decreased operative time
- Improved recovery time
- Works on footswitch or handswitch pencils
- Sterile — individually packaged

Item No.	Description	Qty.
ES-GEAC	Goco Adenoid Curette	1 ea.
ES-GEAC-6	Goco Adenoid Curette	6/bx.

View Adenoid Curette video at www.anthonyproducts.com

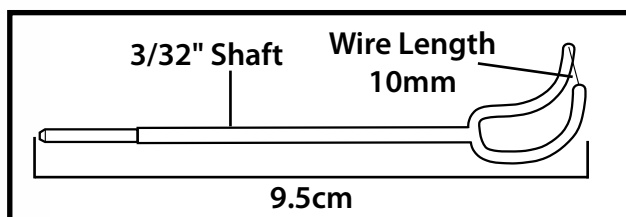

Patent Pending 61/144,841

Goco Silicone Suction Catheter



- Use in conjunction with the Goco Electrocautery Adenoid Curette
- Retracts the soft palate and evacuates smoke during electrosurgical adenoid surgery
- Inserts through the nostrils making the suction catheter completely hands free
- Sterile — individually packaged

Item No.	Description	Qty.
ES-GSSC	Goco Silicone Suction Catheter	1 ea.
ES-GSSC-10	Goco Silicone Suction Catheter	10/bx.

Thermal Cautery Unit

- Provides all the cautery power necessary for treating the numerous minor surgical conditions seen by physicians in everyday practice
- Ideal for dermatological, ophthalmologic, ENT and plastic surgery
- Perfect for any procedure that requires accurate hemostasis
- The autoclavable handpiece, with six-foot cord, accepts an assortment of reusable thermal cautery electrodes, allowing maximum flexibility
- No aberrant radio frequency currents to cause unintended tissue damage, shocking or interference with implantable pacemakers/defibrillators
- Maintains temperature during procedures, unlike battery operated cauteries
- Adjustable knob allows simple and precise levels of cauterization
- TCU unit includes handpiece, ES-G210, ES-G212 and ES-G214 tips

Item No.	Description	Qty.
ES-G150	Thermal Cautery Unit	1 ea. (includes handpiece, A, C & E disposable tips)
ES-G250	Replacement Handpiece with 6ft. cable	1 ea.
ES-G220	Handpiece Sheaths.....	100/bx.
ES-G150-MS	Mobile Stand with Locking Casters	1 ea.

TCU REUSABLE TIPS


ES-G201-A	Loop Tip Cautery	1 ea.
ES-G201-B	Fine Tip Cautery Flat	1 ea.
ES-G201-C	Loop Tip Cautery Angled	1 ea.
ES-G201-D	Ball Point Tip Cautery	1 ea.
ES-G201-E	Fine Point Tip Cautery	1 ea.

TCU DISPOSABLE TIPS

ES-G210	Disposable Thermalite Loop Tip	10/bx.
ES-G212	Disposable Thermalite Loop Tip Angled	10/bx.
ES-G214	Disposable Thermalite Fine Tip	10/bx.



Operates standing upright or laying down

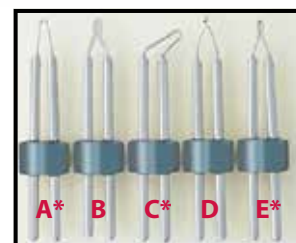

*Tips A, C & E are include with purchase of unit

High & Low Temperature Cauteries

High-Temperature Cauteries

								
AA03	AA05	AA09	AA11	AA17	AA21	AA25	AA27	AA29

High-Temperature Cauteries

- Popular for many types of procedures including pinpoint hemostasis, evacuation of subungual hematoma, controlling diffuse bleeders, vasectomies and sculpting woven grafts
- Individually packaged sterile

Item No.	Description	Qty.	Shelf Life
AA03	High Temp Loop Tip	10/bx.....	4 yr
AA05	High Temp Loop Tip 2" Extended Shaft	10/bx.....	4 yr
AA09	High Temp Loop Tip 5" Extended Shaft	10/bx.....	4 yr
AA11	High Temp Adjustable Fine Tip.....	10/bx.....	4 yr
AA17	High Temp Fine Tip 2" Extended Shaft	10/bx.....	4 yr
AA21	High Temp Vasectomy Tip.....	10/bx.....	4 yr
AA25	High Temp Fine Tip Micro	10/bx.....	2 yr
AA27	High Temp Loop Tip Micro.....	10/bx.....	2 yr
AA29	High Temp Vasectomy Tip Micro.....	10/bx.....	2 yr

Cautery Tip Cleaner

(Item No. XM03-2-00)

- An abrasive foam pad designed to aid in cleaning and removing material from the tips of cautery probes and electro-surgical forceps
- The strong adhesive backside
- 35 per box

WARNING: Not for use with Non-Stick Electrosurgical Forceps!



Low-Temperature Cauteries


			
AA00	AA02	AA04	AA90

Low-Temperature Cauteries

- Specifically designed for ophthalmology
- Fine tips for precise pinpoint hemostasis
- All disposable cauteries are individually packaged sterile

Item No.	Description	Qty.	Shelf Life
AA00	Low Temp Fine Tip	10/bx.....	4 yr
AA02	Low Temp Adjustable Fine Tip	10/bx.....	4 yr
AA04	Low Temp Elongated Fine Tip	10/bx.....	4 yr
AA90	Low Temp Micro Fine Tip	10/bx.....	2 yr

Change-A-Tip Cautery Replacement Tips

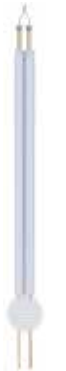										
AA-R100	AA-R101	AA-R103	AA-R104	AA-R105	AA-R106	AA-R109	AA-R112	AA-R110	AA-R111	AA-R121

Item No.	Description	Qty.	Item No.	Description	Qty.
AA-R100	Low Temp Fine Tip.....	10/bx.	AA-R106	High Temp 2" Flex Fine Tip.....	10/bx.
AA-R101	High Temp Fine Tip	10/bx.	AA-R109	High Temp 5" Loop Tip	10/bx.
AA-R103	High Temp Loop Tip.....	10/bx.	AA-R112	High Temp 5" Fine Tip.....	10/bx.
AA-R104	Low Temp Elongated Fine Tip	10/bx.	AA-R110	High Temp Contra Angle Fine Tip	10/bx.
AA-R105	High Temp 2" Flex Tip.....	10/bx.	AA-R111	High Temp Circular Tip.....	10/bx.
			AA-R121	High Temp Vasectomy Tip	10/bx.

Change-A-Tip Deluxe Low-Temp Cautery Kit



Kit includes:

- One Low-Temperature Handle
- Two AA-R100 Non-Sterile Tips
- One "AA" Battery
- One Foam-lined Case

Item No.	Description
AA-RC-L	Change-A-Tip Deluxe Low-Temp Cautery Kit
AA-HIT0	Replacement Handle

Change-A-Tip Deluxe High-Temp Cautery Kit



Kit includes:

- One High-Temperature Handle
- Two AA-R101 Non-Sterile Tips
- Two AA-R103 Non-Sterile Tips
- Two "AA" Battery
- One Foam-lined Case

Item No.	Description
AA-RC	Change-A-Tip Deluxe High-Temp Cautery Kit
AA-HIT1	Replacement Handle

Change-A-Tip Deluxe HI-LO Cautery Kit



Kit includes:

- Two Handles (1 low temp & 1 high temp)
- One Sterile AA-R100 Tip
- One Sterile AA-R101 Tip
- One Sterile AA-R103 Tip
- One Sterile AA-R121 Tip
- Six "AA" Batteries
- One Foam-lined Case

Item No.	Description
AA-DEL2	Change-A-Tip Deluxe HI-LO Cautery Kit

Change-A-Tip Cautery Accessories


Item No.	Description	Qty.
AA-H101-ADH	Fine Tip with Sheath Seal	10/bx.
AA-PT	Protective Tray.....	1 ea.
AA-SS	Cautery Sheaths (sterile)	10/bx.
AA-0007	Subungual Hematoma Nail Drill Bit	1 ea.
	(for use with High Temp Fine Tip Cautery)	



AA-H101-ADH



AA-PT

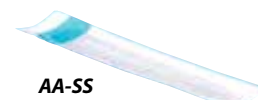

AA-SS



AA-0007

Protect Patients During Electrosurgical Procedures

- Non-conductive insulated coating protects patients and staff from burns and currents
- Blue coating quickly distinguishes from non-coated instruments*
- Autoclavable

**Non-conductive coating on Coated Fraizer Suction Tubes is black*

SPECIAL ORDER INSTRUMENTS

Anthony Products can add a non-conductive coating to most surgical instruments for an additional charge. Please call 800-428-1610 for additional information and pricing.

(Note: special order instruments are non-returnable and cannot be returned for any reason)

Coated Adson Dressing Forceps

- Ideal for removing sutures, holding tissue, moving dressings and removing drainage tubes

Item No.	Description
ES06-20-12SM	4 $\frac{3}{4}$ " , small, smooth jaws



Coated Potts-Smith Forceps

- Designed for grasping fine, delicate tissue
- Narrow tips

Item No.	Description
ES06-46-18	7", straight, serrated
ES06-46-20	8", straight, serrated
ES06-46-24	9 $\frac{1}{2}$ ", straight, serrated, platform
ES06-46-24C	9 $\frac{1}{2}$ ", curved, serrated, platform
ES06-46-25	9 $\frac{1}{2}$ ", straight, serrated
ES06-46-30	12", straight, serrated, platform
ES06-46-80S	8 $\frac{1}{4}$ ", smooth
ES06-46-84S	9 $\frac{1}{2}$ ", straight, smooth
ES06-47-20A	8", angled, platform
ES06-47-24	8", straight, smooth
ES06-47-24A	9 $\frac{1}{2}$ ", angled, platform
ES06-47-25A	10", angled, platform
ES06-47-30	12", straight, smooth



Coated Cushing Forceps

- Used for applying and removing surgical dressings
- Fine, serrated tips

Item No.	Description
ES06-60-17	7", serrated



Coated Gerald Dressing Forceps

- Used for applying and removing surgical dressings
- Fine, serrated tips

Item No.	Description
ES06-60-90	7", serrated



Coated Lucae Bayonet Forceps

- Used in wound dressing and general tissue grasping
- Bayonet style handle offers a less obstructive view

Item No.	Description
ES07-30-14	5 1/2", serrated
ES07-30-19	7 1/2", serrated



Coated Jansen Bayonet Forceps

- Used in wound dressing and general tissue grasping
- Bayonet style handle offers a less obstructive view

Item No.	Description
ES07-30-20	8 1/4", serrated



Coated Micro Jewelers Forceps

- Used for handling fine tissue that can be easily damaged during microsurgical procedures
- Super fine pointed tips
- Flat handles

Item No.	Description
ES07-60-04	4 1/2", straight, size 4
ES07-60-07	5 1/2", curved, size 7



ES07-60-04



ES07-60-07

Coated Crile Double-ended Retractor

- Concave, L-shaped blades
- One small blade and one large blade on opposite end


Item No.	Length	Blade Widths
ES20-09-30	4 1/2"	3 x 11mm/6 x 20mm



Coated Farabeuf Double-ended Retractors

- L-shaped round blades
- Blunt blades prevent tissue damage

Item No.	Description
ES20-10-14	6", small
ES20-10-16	6", large

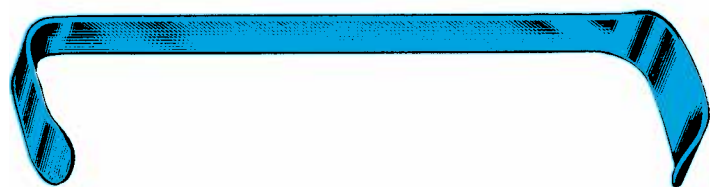

Coated Cottle Hook Retractors

- Used during nasal procedures
- Rounded, ball prongs prevent tissue damage

Item No.	Description	Width
ES20-50-42	6 1/4", double ball prongs	12mm



Insulated Surgical Instruments

Coated U.S. Army Double-ended Retractor

- Commonly used during small wounds and abdominal surgery procedures
- Flat, fenestrated handle
- 90° angled L-shaped blades
- Slightly curved, rounded lips

Item No.	Description
----------	-------------

ES22-18-22	8 1/2", small
------------	---------------

ES22-18-23	8 1/2", large
------------	---------------



Coated Pediatric Deaver Retractors

- Smaller blade widths for use during pediatric surgery
- Used for deep wound retraction during abdominal surgery
- Bottom hook designed to reduce hand fatigue
- Flat, curved blades

Item No.	Length	Blade Width
----------	--------	-------------

ES23-75-00	8 1/2"	5/8"
------------	--------	------

ES23-75-01	7"	3/4"
------------	----	------

ES23-75-02	8"	7/8"
------------	----	------



Coated Adult Deaver Retractors

- Used for deep abdominal wound retraction
- Bottom hook designed to reduce hand fatigue
- Flat, curved blades

Item No.	Length	Blade Width
----------	--------	-------------

ES23-75-03	9"	1"
------------	----	----

ES23-75-04	12"	1"
------------	-----	----

ES23-75-05	13"	1"
------------	-----	----



Coated Deaver Right Angled Retractor

- Used for deep abdominal wound retraction
- Bottom hook designed to reduce hand fatigue
- Right angle blade

Item No.	Length	Blade Width
----------	--------	-------------

ES23-75-23	9"	1"
------------	----	----



Coated Senn Double-ended Retractor

- Used to expose wounds and surface layers of skin
- One end features 90° angled, blunt blade
- Opposite end has a three blunt prong rake tip

Item No.	Length	Width
----------	--------	-------

ES20-05-09	6 1/8"	10mm
------------	--------	------



Coated Desmarres Lid Retractors

- Retracts the upper or lower eye lid
- Keeps eyelashes out of the surgical field
- Curved, saddle style blade


Item No.	Description	Width
ES60-05-11	5 1/8", small	12mm
ES60-05-12	5 1/8", medium	14mm
ES60-05-13	5 1/8", large	16mm
ES60-05-14	5 1/8", x-large	18mm



Coated Fenestrated Desmarres Lid Retractors

- Retracts the upper or lower eye lid
- Fenestrated blade allows an unobstructed view of the operating site
- Curved, saddle style blade

Item No.	Description	Width
ES60-05-22	5 1/8", medium	14mm
ES60-05-24	5 1/8", x-large	18mm

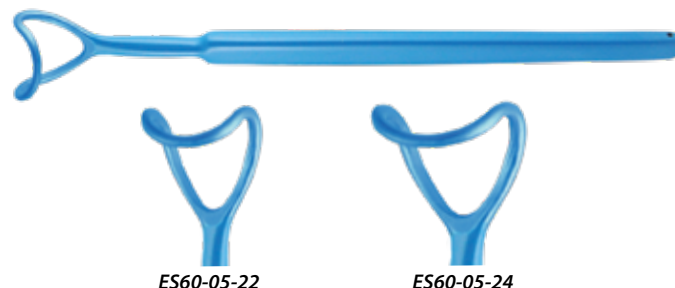

Coated Jaeger Lid Plate

- Protects the cornea and globe during eyelid procedures
- One end is concave and other end is convex
- One blade size 24mm, other side is 20mm

Item No.	Description
ES60-05-95	4", lid plate



Coated Weider Tongue Depressors

- Used for retracting and protecting the tongue during oral surgical procedures
- Serrated blade with fenestrations allow fluids to be removed with suction tubes
- L-shaped curved handle
- Bottom hook to reduce hand fatigue

Item No.	Description	Blade Width
ES72-01-40	5 1/2", small, serrated	28mm
ES72-01-41	5 1/2", large, serrated	30mm



Coated Fraizer Suction Tubes

- Used to remove fluid and object from ear or nasal passages
- Thumb plate with finger cutoff for controlling suction
- Includes style to remove debris inside of suction tube

Item No.	Length	Tip Diameter
ES09-83-08	4"	8 French (2.67mm)
ES09-83-10	4"	10 French (3.33mm)
ES09-83-12	4"	12 French (4mm)
ES09-83-08-L	6"	8 French (2.67mm)
ES09-83-10-L	6"	10 French (3.33mm)
ES09-83-12-L	6"	12 French (4mm)



**Frazier Suction Tubes
require this cable**
(Item No. ESFS-C)
Length: 10ft.



ELECTROSURGICAL EQUIPMENT

P R O D U C T C A T A L O G

Anthony Products, Inc.

7740 Records Street • Indianapolis, IN 46226

(317) 545-6196 • 800-428-1610 • Fax (317) 543-3289

www.anthonyproducts.com

LNSEBFVER08-03-11-21